

**With secured results against COVID-19**

**ViroQ SARS-CoV-2 - Real-time PCR Kit**

# Your benefits

## One-step multiplex real-time SARS-CoV-2 Test

The ViroQ SARS-CoV-2 Kit is designed for specific detection of viral RNA Gene E and RdRP of SARS-CoV-2 in respiratory samples

## Be confident of negative results

- False negative results due to a lack of cell material in the sample are excluded by simultaneous detection of a universally expressed human gene
- A mix of human RNA and SARS-CoV-2 reference RNA is used as positive control (included in the kit)

## Our product features - Your benefits

- Highly specific and exclusive detection of the SARS-CoV-2 genes E and RdRP
- One-Step real-time PCR assay - easy to set up
- Results in 90 minutes – save time
- Works with common RNA-Isolation kits – be independent from single suppliers
- Ambient shipping condition - no extra costs and handling of dry ice
- Quality made in Germany – high supply chain security

## Compatible with most real-time cyclers


- The multiplex real-time PCR kit uses three colours that are compatible with most common real-time cyclers on the market
- Tested on QuantStudio 6 (Thermo Fischer Scientific), CFX 96™ (Bio-Rad) Lightcycler® 480 (Roche)

# Workflow and Interpretation



## Example for positive result

**RdRP gene:** FAM  
**E gene:** CFO560 / HEX / VIC / JOE  
**Internal control:** CFR610 / Texas Red / ROX

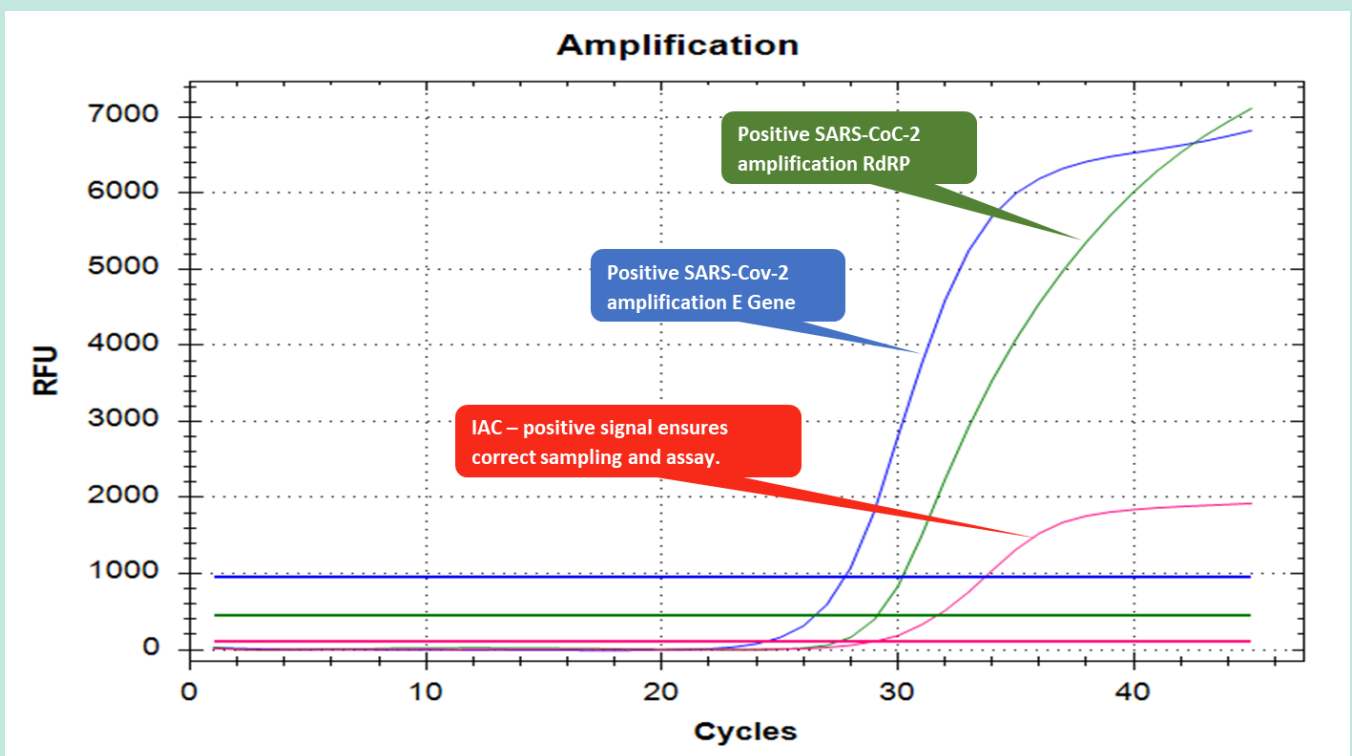

## Ordering Information

Cat. No.	Kit	Size
728250	ViroQ SARS-CoV-2	96 Tests
728251	ViroQ SARS-CoV-2	480 Tests

**BAG Diagnostics GmbH**  
Amtsgerichtsstr. 1-5  
35423 Lich / Germany

### Ordering:

Phone: +49 (0) 6404 925 450  
Fax: +49 (0) 6404 925 460  
Email: [order@bag-diagnostics.com](mailto:order@bag-diagnostics.com)

### Customer Service:

Phone: +49 (0) 6404 925 125  
Fax: +49 (0) 6404 925 421  
Email: [info@bag-diagnostics.com](mailto:info@bag-diagnostics.com)