

**Easily increase test throughput:
In one test: COVID-19 Flu A/B RSV**

ViroQ SARS-FluA/B-RSV RT-PCR Kit

**COMING
in
October**

Your benefits

One-step multiplex real time SARS-CoV-2, Influenza and RSV Test

The ViroQ SARS-FluA/B-RSV - Kit is designed for specific detection of viral RNA genes N and RdRP of SARS-CoV-2, Influenza A and separate B as well as RSV in respiratory samples.

Be confident of negative results

- False negative results due to a lack of cell material in the sample are excluded by simultaneous detection of a universally expressed human gene
- A mix of human RNA and SARS-CoV-2, Influenza and RSV reference RNA is used as positive control (included in the kit)

Our product features - Your benefits

- Highly specific and exclusive detection of the SARS-CoV-2 RdRP and N gene, Influenza A, Influenza B and respiratory syncytial virus
- One-Step real time PCR assay - easy to set up
- Results in approximately 60 minutes – save time
- Works with common RNA-Isolation kits – be independent from single suppliers
- Ambient shipping condition - no extra costs and handling of dry ice
- Quality made in Germany – high supply chain security

Compatible with most real time cyclers


- The multiplex real time PCR kit uses five colours. Four for the specific genes and a cycler-specific colour for the detection of the universally expressed human gene, so that the test is compatible with most common real time cyclers on the market.
- Validated on QuantStudio 6 (Thermo Fischer Scientific), CFX 96™ (Bio-Rad) Lightcycler® 480 (Roche)

Workflow and Interpretation



Example for positive result

RdRP and N gene; Flu A; Flu B; RSV
Included Separate Internal control for: CFX / LightCycler® 480 / QuantStudio™

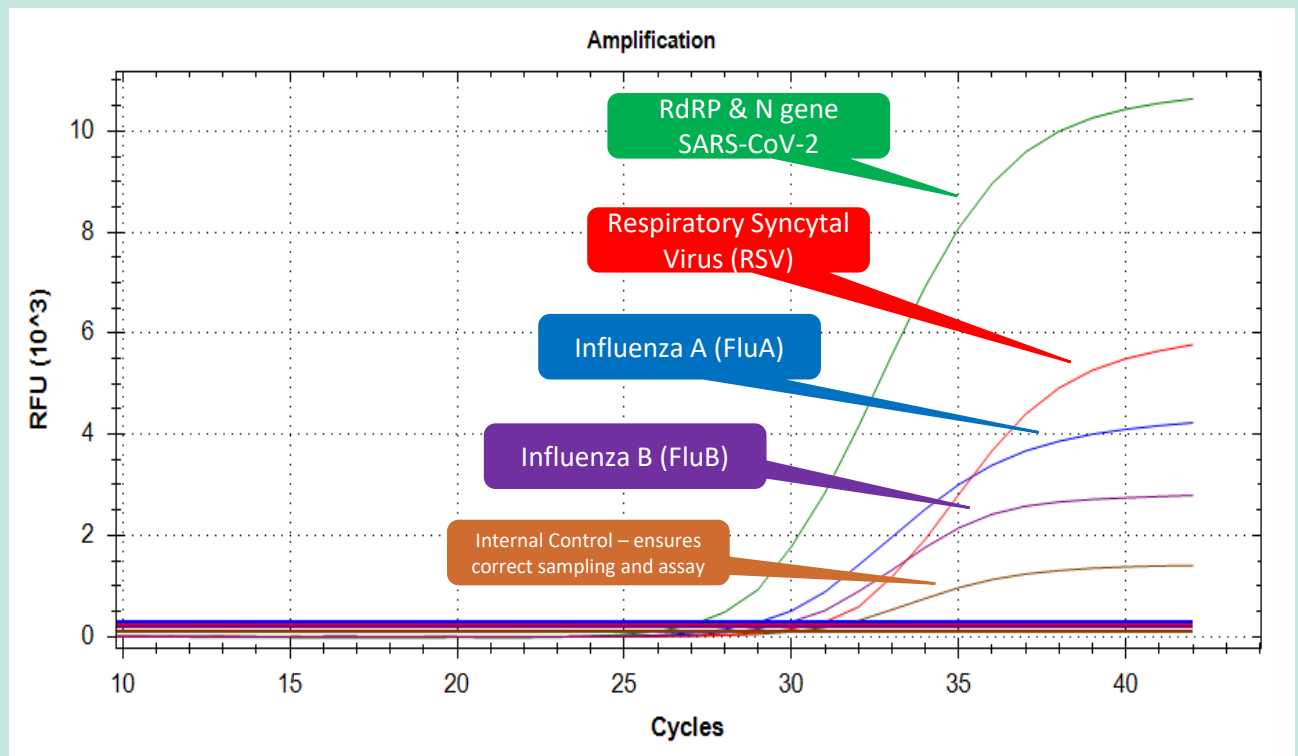

Artificial sample to show all Signals; this case has an extremely low probability of occurrence

Ordering Information

Combination test – for SARS-CoV-2, FluA, FluB, RSV

Cat. No.	Kit	Size
728267	ViroQ SARS-FluA/B-RSV	96 tests
728268	ViroQ SARS-FluA/B-RSV	480 tests

COMING
in
October

SARS-CoV-2 test

Cat. No.	Kit	Size
728250	ViroQ SARS-CoV-2	96 tests
728251	ViroQ SARS-CoV-2	480 tests
728263	ViroQ Rapid SARS-CoV-2	96 tests
728264	ViroQ Rapid SARS-CoV-2	480 tests

Ordering:

Phone: +49 (0) 6404 925 450
Fax: +49 (0) 6404 925 460
Email: order@bag-diagnostics.com

Customer Service:

Phone: +49 (0) 6404 925 125
Fax: +49 (0) 6404 925 421
Email: info@bag-diagnostics.com